




COVER-DRI


Cover-Dri is the Attends range of surface protection underpads, specially designed to manage light to moderate urinary incontinence and for clinical procedures where non-sterile containment is required.




Fitting Guide

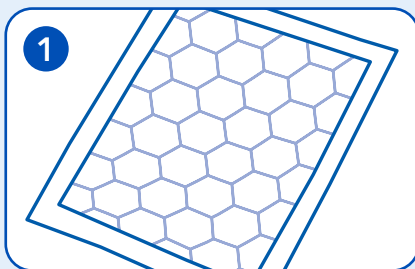


Fully open the underpad ensuring the waterproof backsheet is positioned on the surface you are positioning the person and the soft white absorbent cushions are next to the skin.

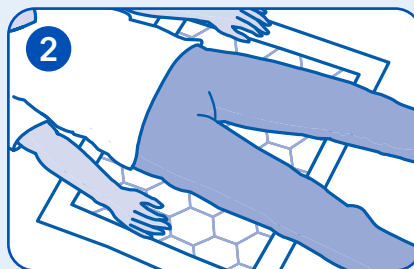


Position the person on the underpad.

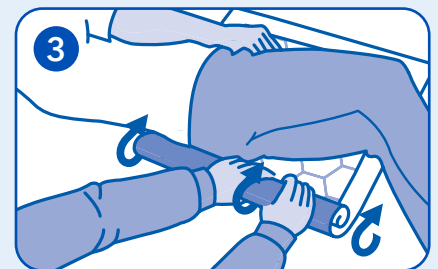


Turn the person into the lateral position and roll up the used underpad towards the person's back.

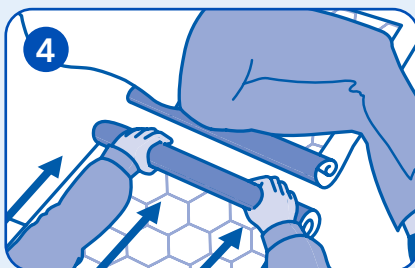

Unfold the new underpad, positioning it towards the back of the person.

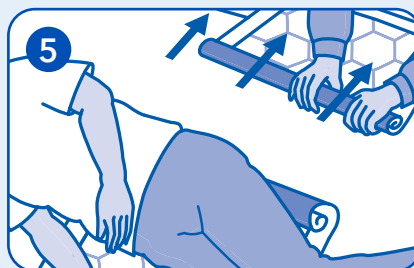

Turn the person on to the underpad, ensuring the pad is lying flat. Remove the used product.

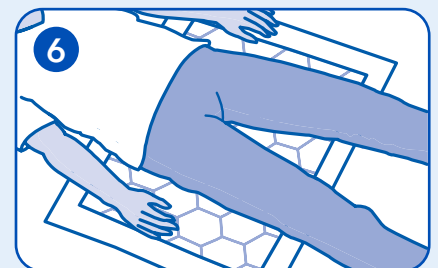

Lay the person in their original position.

COVER-DRI

Hints and Tips

Selection

- Choose the correct product size and absorbency depending on the clinical situation you are using the product in.

Disposal

- Dispose as per your local authority's waste policy.
- Do not flush the pad down the toilet.

Leakage

- Check you have the correct size and absorbency of underpad for the situation it is being used in.



Skin tolerability is dermatologically approved by SGS proderm Institute for Applied Dermatological Research.



Attends Healthcare AB
Järnvägsgatan 4
SE-57833 Aneby
Sweden



At Attends, we care about each other,
our community and the planet.



Products produced by using renewable energy and supporting verified emission reduction projects.



These products have achieved the certification of the recognised Nordic Swan Ecolabel.



Plastic bags made with min. 30% recycled material.



The cellulose within these products originates from sustainably managed forests and controlled sources.



With this label these products are certified as being harmless for human health.