




SOFT


Attends Soft is our range of discreet, absorbent liners and pads that are kind to the skin and can be worn in normal underwear. With reliable protection for very light to heavy urinary incontinence, they help restore confidence, comfort and a sense of normality in everyday life.




Fitting Guide

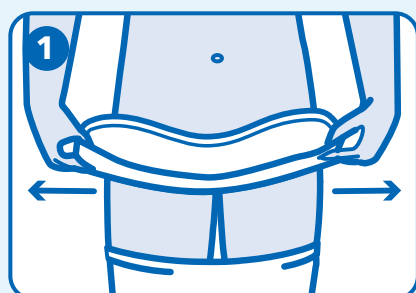


Open the product fully, then cup and fold it to activate the leakage barriers. This does not apply to Soft 0, 1 and 1 Long.

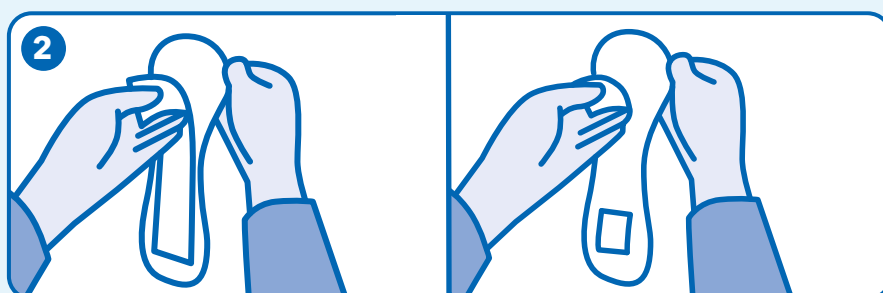

Depending on the product's absorbency, there will be one or two fixation tapes (0 - 4 has one, 5 - 7 has 2). To fit the product, remove the adhesive tape on the back of the pad.

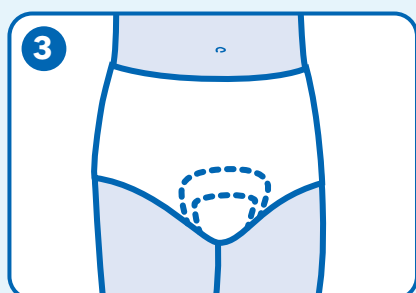

Place the pad securely in your own close fitting underwear.

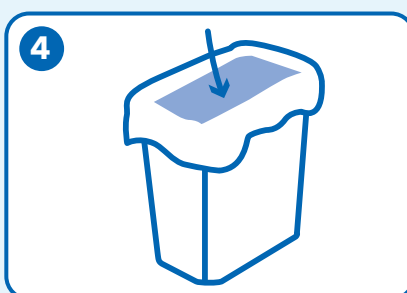

After use, roll up the used pads and dispose of them in general waste.

SOFT

Hints and Tips

Selection

- Please select the correct product absorbency level.
- For optimal performance, we suggest using the product in combination with close fitting underwear or with Attends Stretch Pants for a closer fit.
- The smaller the pad chosen, the greater the discretion and comfort for the user, whilst also helping to reduce waste.

Disposal

- Roll up pads and dispose of as per your local authority or NHS waste policy.
- You may use the single wrap, if provided, to dispose of the product.
- Do not flush the pad down the toilet.



Skin tolerability is dermatologically approved by proDERM Institute for Applied Dermatological Research